manual


Massage gun

Art no: 47900

EN

NO

SV

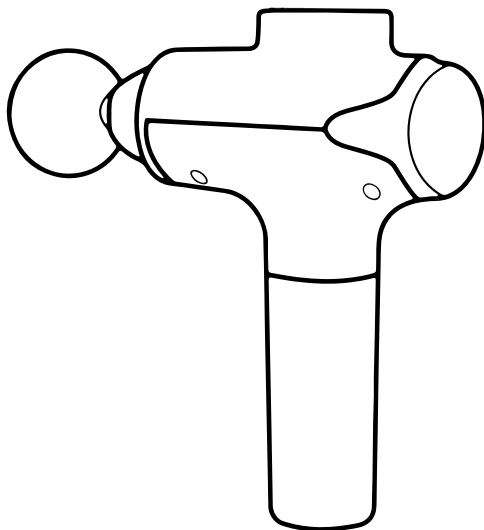www.rubicson.com
Box 50435,
Malmö, Sweden
2024-07-05


Rubicson®


Overview

Message gun

1. Massage head piston
2. Battery release button
3. Battery/grip
4. LED indicator
5. Speed level indicator
6. Decrease speed level
7. Battery indicator
8. Increase speed level
9. Charging port
10. Power on/off

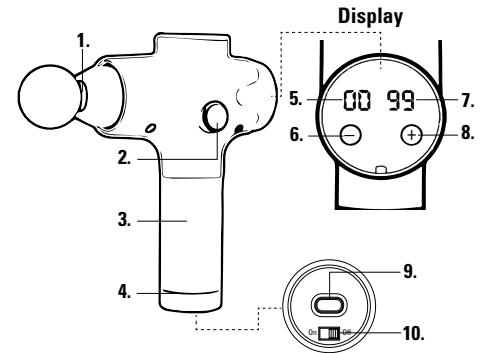Message heads

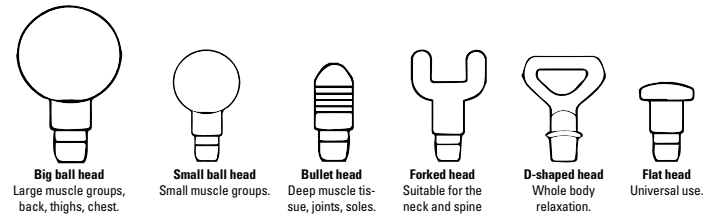Specifications

Motor:	20 W brushless
Turn speed:	1200–2700 rpm
Amplitude:	12 mm
Speed settings:	20
Battery life:	3–4 hrs
Battery:	11.1 V lithium-ion (1,500 mAh)
Charge time:	Approx. 4 hrs
Noise level:	50 dB (max.)
Included:	Massage gun, 6 massage heads, USB charging cable, case, manual

Use

Charging

Connect the included USB charging cable between the charging port (9) on the massage gun and a USB charger (sold separately). The LED indicator (4) lights up and flashes red during charging, then changes to flashing green when the battery is almost fully charged (>90% battery capacity), and stays lit green when the battery is fully charged. The battery indicator (7) shows the current battery level in real-time during use (from 20–99 %).

Note! Do not use the massage gun during charging. Disconnect the charging cable when the battery is fully charged. We recommend turning the product off during charging.

Mounting a massage head

Attach the massage head you want to use to the massage gun piston (1) by pressing the bottom of the head into the front of the massage gun. Make sure the head is securely fastened before moving to the next step.

Power on/off and speed settings

Set the Power switch (10) to "On" to power on the massage gun. The LED indicator lights up and shows the remaining battery capacity: green at 99–20% and red when the capacity is below 20%. The display shows the current battery level (7) and massage speed setting (5) (preset to "00"). Press the head of the massage gun to the muscles you want to massage.

Press the buttons (6, 8) to increase and decrease the speed level to adjust the speed. It is recommended to start with a low setting and slowly increase the speed for a more enjoyable massage. There are 20 speed levels available where level 1 runs at 1200 rpm and level 20 runs at 2700 rpm.

Press and hold both speed buttons (6, 8) for 3 seconds to stop the massage gun and reset the speed level to 0.

Set the Power switch (10) to "Off" to turn off the massage gun and place it back in the case before storing it.

Note! The massage gun automatically stops after 10 minutes of continuous use or if too much force is used during a massage. Press the button to increase speed level (8) to restart the massage. To ensure the longevity of the motor, turn the machine off after 1 hour of use and let it rest for 30 minutes before using the massage gun again.

To increase the battery life further, don't forget to turn the massage gun off when you've finished the massage or during charging.

Maintenance

Cleaning

Clean the massage gun and the massage heads using a damp towel and dry them with a soft cloth. Do not use any detergents on any of the parts.

Battery

Note! The massage gun battery is a lithium-ion battery. Over time and with regular use, the chemical composition and high energy output of any lithium-ion battery naturally diminish, resulting in a decrease in performance and capacity. A new battery can be purchased and replace the original battery if the capacity is too low, contact the place of purchase if you want to buy a new battery.

Remove the battery from the massage gun if you are not going to use it for long periods of time. Make sure the Power on/off switch (10) is set to "Off". Press the Battery release button (2) on the side of the massage gun and pull the battery grip down to separate it from the massage gun.

Make sure you align the battery correctly before reinstalling it as it only fits in a single position. Align the bulge on the battery with the battery release button and press the battery into place.

Safety information

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

Children shall not play with the appliance. User maintenance shall not be made by children without supervision.

Do not use the massage gun on injured body parts before consulting a doctor.

Never attempt to disassemble the product. Repairs should only be performed by a qualified technician. By repairing it yourself there is a risk of electric shock and/or short circuit.

Keep the device away from moisture, water, high temperatures, and direct and strong sunlight.

Do not use the massage gun to directly massage dry skin to avoid injuries. The massage area should either be separated from the massage gun by a dry towel or clothing or used with massage oil. It is not recommended to massage a single muscle or part of the body for more than 60 seconds at a time.

Only use high speeds and high-pressure settings on soft tissue parts of the body to avoid causing pain and discomfort. Do not use high speed and high pressure on the head and any bone parts of the body. Frequent use of high speeds and high pressure on the same body part may cause abrasions. If you feel pain and/or discomfort, stop using the massage gun immediately.

Keep fingers, hair and/or any other part of the body away from the moving telescopic rod of the massage gun and the heat dissipation vent of the motor to prevent injuries.

Do not cover the heat dissipation vent or insert any foreign objects into the motor.


If the power adapter is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified personnel to avoid hazard. Only use the included power adapter.

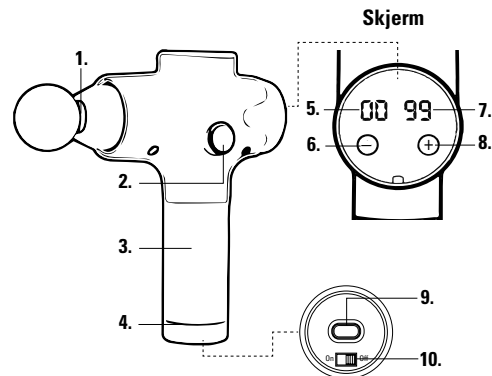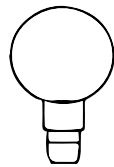
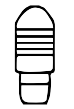
Do not leave the massage gun running unsupervised.

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