manual


Daylight therapy lamp

Item: 64804

EN

NO

SV

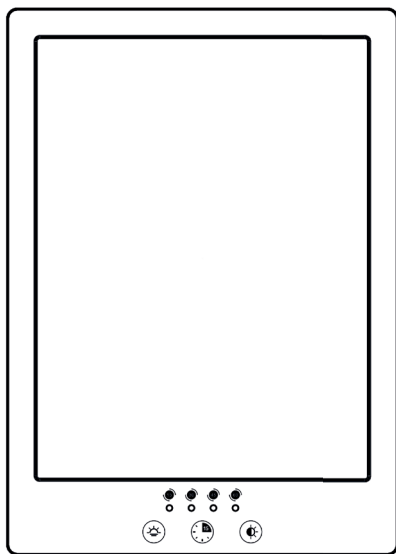
www.cleverio.se
Box 50435
Malmö, Sweden
2021-0714

Cleverio®

Overview

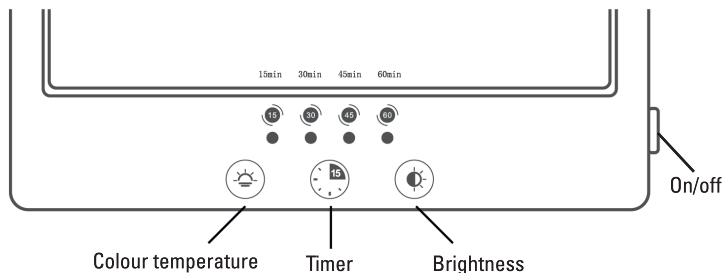Specifications

Power supply	20 V, 0.9 A
Power	3-18 W
Colour temperature	3000-6500 K
Brightness	2000-12000 lux
Timer	15/30/45/60 min
Size	249x178x10 mm
Weight	460 g

Use

Connect the power adapter to the power port on the back of the lamp and to an outlet. Press the On/Off button to turn on the lamp.

Place the lamp 15-60 cm from you. Place it in a 45-degree angle in front of you.

Press and hold the Brightness button until the lamp reaches the desired brightness. Consider using a brighter light in the morning, and a darker light in the afternoon.

Press and hold the Colour temperature button until the lamp reaches the desired colour temperature.

Press the Timer button to set the timer. Press repeatedly to set the timer to 15, 30, 45 or 60 minutes. The corresponding indicator lights up. If no timer is set, the lamp shines continuously.

Recommended use

You do not need to sit in front of the lamp for a long time but rather use it for 30-60 minutes every day or split it up into several shorter sessions during the day. Try to use the lamp regularly and at the same time every day for maximum effect.

Safety information

Avoid using the lamp near water.

The lamp is for indoor use only.

Children and others who may not understand the need for these precautions may only use the lamp under close supervision.

Certain medications (e.g. anti-depressants) and diagnosed mental health conditions may induce light sensitivity, which may result in discomfort, hyperactivity, or hypomania. Consult with your health care provider before using the lamp.

Consult with your health care provider if you have a pre-existing eye condition or are at risk of developing age-related macular degeneration or other eye related conditions.

Stop using the lamp if you experience any discomfort.


This product is not to be used with any kind of light dimmers, timers, motion detectors, voltage transformers or extension cords.

Do not use the product if the lamp or power cord has been damaged.

Do not cover the lamp during use.

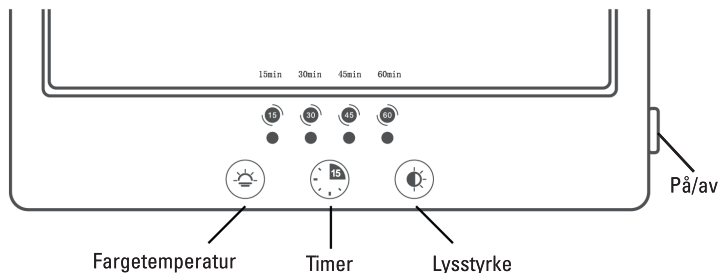Do not clean the lamp with abrasive solvents or cleaners. Only clean the lamp with a dry, soft cloth.

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